

PRODUCT INFORMATION

Common Name TJ-011133, TJC-4 Conjugate Unconjugated CD47;MER6;IAP **Synonyms Applications** ELISA; Flow Cyt

Recommended

Dilutions

ELISA 1:5000-10000; Flow Cyt 1:100

Formulation & Reconstitution Lyophilized from sterile PBS, pH 7.4. Normally 5 % – 8% trehalose is added as protectants before lyophilization. Please see Certificate of Analysis for specific instructions of reconstitution.

Host Species Humanized

IgG type IgG4 Reactivity Human **Target CD47 Uniprot ID** Q08722

Anti-CD47(lemzoparlimab biosimilar) mAb **Description**

Delivery In Stock

> Store at -20°C to -80°C for 12 months in lyophilized form. After reconstitution, if not

intended for use within a month, aliquot and store Storage & Shipping

at -80°C (Avoid repeated freezing and thawing). Lyophilized proteins are shipped at ambient

témperature.

Research grade biosimilar. Not for use in

Background therapeutic or diagnostic procedures for humans

or animals.

Usage Research use only

> All DIMA recombinant antibodies are genuinely generated by DIMA Biotech. They are all under patent application. Any protein sequencing or reverse engineering attempt is prohibited. We are

DIMA Disclaimer

actively scrutinizing all patent application to ensure no IP infringement.

Email: info@dimabio.com Website: www.dimabio.com



Anti-CD47 (lemzoparlimab biosimilar) mAb ELISA

0.1 μg of Human CD47, mFc-His tagged protein per well

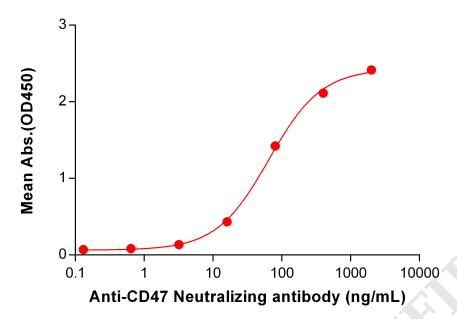


Figure 1. ELISA plate pre-coated by 1 μ g/mL (100 μ L/well) Human CD47 protein, mFc-His Tag PME100008 can bind Anti-CD47 Neutralizing antibody (BME100076) in a linear range of 3.2-400 ng/mL.

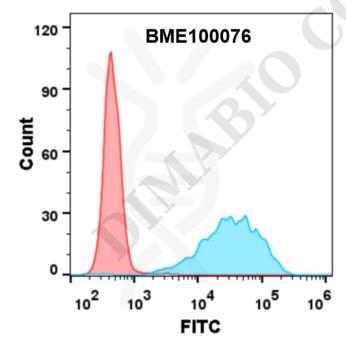


Figure 2. Jurkat cell line were surface stained with anti-CD47 (BME100076) (Blue histogram) and isotype control (Red histogram) 1µg/ml followed by Alexa 488-conjugated Goat anti-Human IgG secondary antibody.

Email: info@dimabio.com Website: www.dimabio.com





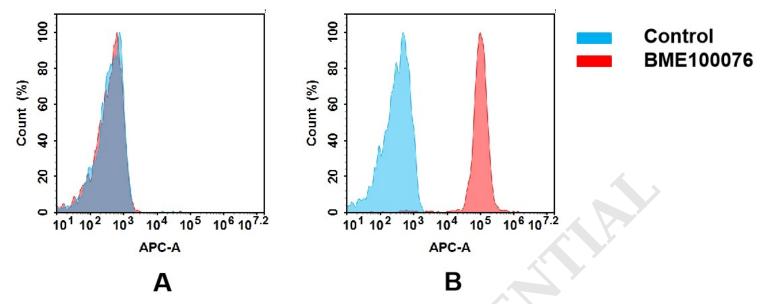


Figure 3. Flow cytometry analysis of antigen binding of anti-human CD47 mAb(BME100076). (A) BME100076 does not bind to CHO-S cells that do not express CD47. (B) A clear peak shift of BME100076 was seen compared to the control when incubated with CD47-expressing 8226 cells, indicating strong binding of BME100076 to CD47. Antibodies were incubated at 5 μ g/mL.



