

PRODUCT INFORMATION

Tag	N-MBP Tag, C-Flag&Strep Tag
Target	TRPV1
Synonyms	VR1
Description	Human MBP-TRPV1-Strep full length protein-synthetic nanodisc
Delivery	In Stock
Uniprot ID	Q8NER1
Expression Host	HEK293
Protein Families	Druggable Genome, Ion Channels: Transient receptor potential, Transmembrane
Protein Pathways	Neuroactive ligand-receptor interaction
Molecular Weight	The human full length MBP-TRPV1-Strep protein has a MW of 135.0 kDa Lyophilized from nanodisc solubilization buffer (20 mM Tris-HCl, 150 mM NaCl, pH 8.0). Normally 5% - 8% trehalose is added as protectants before lyophilization. Please see Certificate of Analysis for specific instructions. Lyophilized from PBS.
Formulation & Reconstitution	Normally 5% - 8% trehalose is added as protectants before lyophilization. Please see Certificate of Analysis for specific instructions. Do not use solvents with pH lower than 6.5 in subsequent experiments.
Storage & Shipping	Store at -20°C to -80°C for 12 months in lyophilized form. After reconstitution, if not intended for use within a month, aliquot and store at -80°C (Avoid repeated freezing and thawing). Lyophilized proteins are shipped at ambient temperature.
Background	Capsaicin, the main pungent ingredient in hot chili peppers, elicits a sensation of burning pain by selectively activating sensory neurons that convey information about noxious stimuli to the central nervous system. The protein is a receptor for capsaicin and is a non-selective cation channel that is structurally related to members of the TRP family of ion channels. This receptor is also activated by increases in temperature in the noxious range, suggesting that it functions as a transducer of painful thermal stimuli in vivo. Four transcript variants encoding the same protein, but with different 5' UTR sequence, have been described for this gene.
Usage	Research use only
Conjugate	Unconjugated



ELISA assay to evaluate MBP-TRPV1-Strep-Nanodisc 0.2 μ g Human MBP-TRPV1-Strep-Nanodisc per well

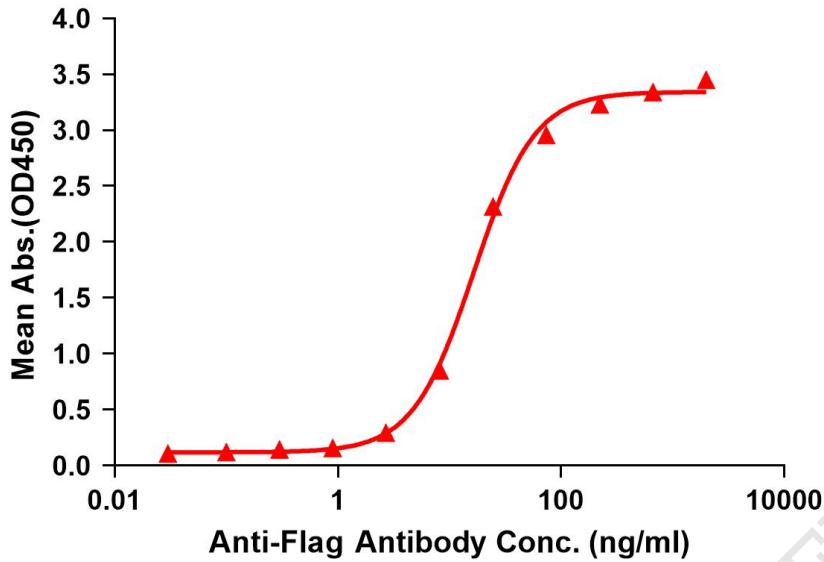


Figure 1. Elisa plates were pre-coated with N-MBP Tag, C-Flag&Strep Tag MBP-TRPV1-Strep-Nanodisc (0.2 μ g/per well). Serial diluted anti-Flag monoclonal antibody solutions were added, washed, and incubated with secondary antibody before Elisa reading. From above data, the EC50 for anti-Flag monoclonal antibody binding with MBP-TRPV1-Strep-Nanodisc is 16.57ng/ml.

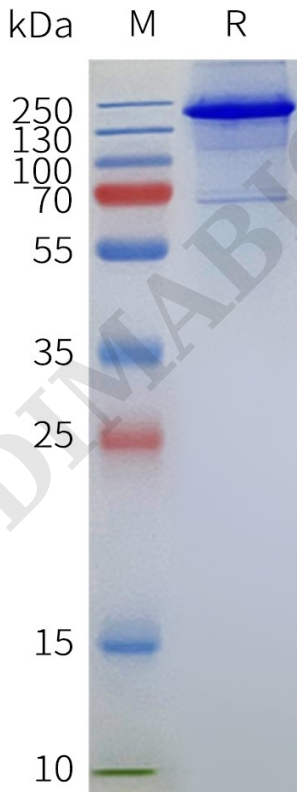


Figure 2. Human MBP-TRPV1-Strep-Nanodisc with N-MBP Tag, C-Flag&Strep Tag on SDS-PAGE

